





The **Invacare® Solara® 3G Tilt-in-Space Wheelchair** offers more adjustability and is easier to tilt as compared to the previous generation Solara wheelchair. With a 300 lb. standard weight capacity, tilt assist adjustment and a seat pan that grows, the Solara 3G wheelchair offers everything you loved about the Solara 2G wheelchair plus much more.

Features

The unique Tilt Assist feature makes tilting easier by reducing the effort required to tilt an occupied wheelchair With three settings, front, middle and rear, the simple Tilt Assist adjustments help to correct for uneven weight distribution

The innovative Low Shear Recline minimizes shearing forces and prevents displacement of positioning tools such as laterals and headrest. Used in combination with tilt, recline can offer additional functional and clinical benefits to users and caregivers

Foot-Operated Wheel Lock easily locks both rear wheels by pressing 1 of the 2 levers of the wheel lock

20" and 24" back canes include a 10° bend option accommodating after-market back hardware while the back angle adjustment location better fits after-market seating systems. In addition, adjustable height back canes can be ordered with or without adjustable angle stroller handles

The DeGage Dynamic Rocker Back accommodates those with involuntary thrusts by helping to reduce pain or discomfort for the user. See www.degage.us for complete details

For more detailed information visit www.invacare.com/solara3g

HCPCS Code: E1161












- 1. Low shear recline
- 2. Back cane options
- 3. DeGage rocker back
- 4. Foot-operated wheel lock option
- 5. Transport Ready Option (TRRO)*
- 6. Transport Brackets (TRBKTS)



*TRRO includes four factory-installed brackets and wheelchair anchored pelvic belt. TRRO has been crash-tested in accordance with ANSI/RESNA WC Vol 1. Section 19 Frontal Impact Test requirements for wheelchairs with a 168 lb. crash dummy, which corresponds to a person with a weight of 114 lb. to 209 lb. TRBKTS includes four factory-installed wheelchair transport brackets. TRBKTS has not been crash tested in accordance with WC19. Use these transport brackets only to secure an unoccupied wheelchair during transport. As of this date, the Department of Transportation has not approved any tie-down systems for transportation of a user while in a wheelchair, in a moving vehicle of any type. It is Invacare's position that users of wheelchairs should be transferred into appropriate seating in vehicles for transportation and use be made of the restraints available by the auto industry. Invacare cannot and does not recommend any wheelchair transportation system

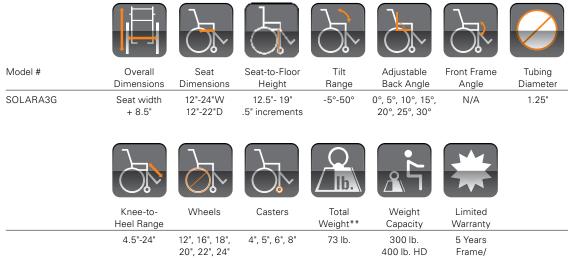

Invacare Corporation www.invacare.com

USA

One Invacare Way Elyria, Ohio 44035-4190 (800) 333-6900

Canada

570 Matheson Blvd. E., Unit 8 Mississauga, Ontario L4Z 4G4 Canada (800) 668-5324

©2012

Invacare Corporation. All rights reserved. Trademarks are identified by the symbols [™] and B. All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted. Specifications are subject to change without notification. Form No. 12-248 Rev. 09/12 120363



** 18x18 frame with 12" wheels

Crossbraces Please see the Owner's Manual and Invacare Configurator for complete option/specification combinations. All references to HCPCS codes: Providers are responsible for determining the appropriate billing codes when submitting claims to the Medicare program, and should consult an attorney or other advisor to discuss specific situations in further detail