



Molift Quick Raiser 205

Ergonomic sit-to-stand lift with
outstanding transfer capability

 **etac**[®]

Molift Quick Raiser 205

The Molift Quick Raiser 205 is a state of the art sit-to-stand lift with outstanding transfer capabilities and exceptional maneuverability. The design encourages a natural movement pattern and fulfils the lifting needs of users with some weight-bearing capacity.

The natural movement pattern

The inclined column encourages a natural movement pattern and fulfils the lifting needs of users up to 452 lbs, with some weight-bearing capacity. The user can feel safe and comfortable in every movement as the lifting experience is the same regardless of the carer, lifting arm or sling.

Exceptional maneuverability

The Quick Raiser 205 is compact, stable and exceptionally easy to maneuver. The lift has a small turning radius and good maneuverability even in tight spaces, thanks to the unique wheel position, angled push rod and outstanding weight distribution of the base.

The lift fits under most beds and chairs, and the variant with electrically adjustable legs allows accessibility almost anywhere. The fixed leg option enables a very small footprint.



Flexible lifting arm

To facilitate the carer's work and not worry stressed users, the lifting arm can be folded away from the user towards the column when not in use.



Excellent accessibility

Quick Raiser 205 fits under most beds and chairs, and the electrically adjustable legs enable accessibility almost everywhere.



Easily maneuvered

Easy to maneuver even in tight spaces thanks to the base's small turning radius and weight distribution.



Ergonomic for the carer

The adjustable legs can remain closed, regardless of movement, giving the carer room to work with minimized risk of injury.

Ergonomic transfers with a natural movement pattern

Consistent & safe lifting

The lifting trajectory follows a natural movement pattern due to the design of the inclined column.

Thanks to this natural movement pattern, the sling stays in place and does not slide up, which is a usual problem with sit-to-stand lifts. As a result, the lifting path is always the same regardless of the sling and who the carer is. This results in a safe, comfortable, consistent lifting experience for the user.



Watch the video for a product presentation.

Molift QuickRaiser 205

Model item no.

Variant	Battery	Castors	Adjustable leg	Fixed leg (FL)
Active, 4-point lifting arm	SLA	Small	-	M29150
Active, 4-point lifting arm	SLA	Large	M29100	-
Active, 4-point lifting arm*	SLA	Large	M29100-00	M29200-00

*alternative leg support cushions (large) with leg straps

For complete variant matrix, contact your local distributor.

Dimensions

Height	44.9 inch
Length (FL)	52 inch (45.9 inch)
Chassis external width	24.8 inch
Legs closed internal/external width	8.7/13.8 inch
Legs open internal width	23.2 inch
Footboard height:	
Castor size small	2.6 inch
Castor size large	3.5 inch

Weight

Maximum user weight SWL	452 lbs
Total weight (FL)	98 (93.7 lbs)

Technical

Battery	SLA 24V 2,9 Ah
Charge time	5 hours
Rated performance	100 lifts/165 lbs
Turning radius with closed/open legs	45 inch/51 inch
Hoisting speed	1.4 inch/sec/165 lbs



Clean / disinfect the product using a standard non-abrasive cleaning agent with a pH between 5 and 9, or with a 70% disinfectant solution.



Accessories & customization

With a wide range of accessories the Quick Raiser 205 is highly customizable

Alternative lifting arms



Active lifting arm
For RgoSling Active slings.



Loop lifting arm
For loop slings




StandUp lifting arm
For RgoSling Standup slings.



Clip sling lifting arm
For clip slings

Slings

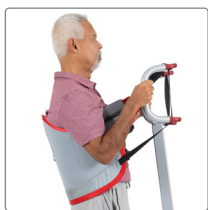

RgoSling Active sling
The 4-point sling is available in sizes XS to XL.

For item numbers and measurements see www.etac.com



RgoSling StandUp sling
The 2-point sling is available in sizes XS to XL. Additional variant is available when more pelvic support is needed.

For item numbers and measurements see www.etac.com

Leg supports/cushions



Bracket & Cushions kit
The alternative leg supports benefits an option of small or large cushions, and a leg strap for safety measurements

Kit with bracket & small cushions ... 2910059
Kit with bracket & large cushions 2910056



Set of Cushions
Only contains cushions, no leg straps.

Set of small cushions..... 2910076
Set of large cushions 2910077

Footboard/inlay



Footboard
This accessory benefits a 10 cm (4inch) elevated footrest. The attachment key is easy accessible and makes it convenient to lock or unlock the accessory to the Quick Raiser.

Footplate..... 2910065



Set of inlay
Inlays for footplate and footboard with slip-resistant material to increase friction.

Set of inlay 2920151

For the latest news and continuously updated product information – visit: www.etac.us.com