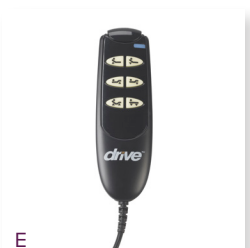
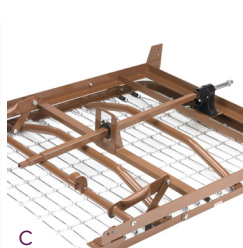


Delta™ Ultra- Light 1000, Semi-Electric Bed

Our Most Popular Model Became
Even Better...**The Game Changer**

**accessories****15030HBC**
Sold Separately**15030HBL**
Sold Separately

Item #	Product Description	UOM
15030	Ultra-Light Semi-Electric Bed	1/cs
15030BV-FR	Ultra-Light Semi-Electric Bed with Full-Length Side Rails (item #15001ABV)	1/cs
15030BV-HR	Ultra-Light Semi-Electric Bed with Half Rails (item #15201BV)	1/cs
15030BV-PKG	Ultra-Light Semi-Electric Bed with Full-Length Side Rails (item #15001ABV) and 80" Innerspring Mattress (item #15006)	1/cs
15030BV-PKG-T	Ultra-Light Semi-Electric Bed with Full-Length Side Rails (item #15001ABV) and 80" Therapeutic Support Mattress (item #15019)	1/cs
15030BV-PKG-1	Ultra-Light Semi-Electric Bed with Half Rails (item #15201BV) and 80" Innerspring Mattress (item #15006)	1/cs
15030BV-PKG-1-T	Ultra-Light Semi-Electric Bed with Half Rails (item #15201BV) and 80" Therapeutic Support Mattress (item #15019)	1/cs
15030BV-PKG-2	Ultra-Light Semi-Electric Bed with Full-Length Side Rails (item #15001ABV) and 80" Foam Mattress (item #15007)	1/cs
15096-SEMI	Ultra-Light Semi-Electric Bed Reclamation Kit	1/ea

- The Delta 1000 is truly universal because the bed frame can be used with Drive's and most other manufacturers' bed ends, old and new.
- The transition box mounted on the foot section is key to interchangeability. (Figure A)
- The transition box allows you to change the rotation of the high/low shaft.
- Headboard is taller than footboard to better fit home decor. (Figure B)
- The wood grain panels are break and scratch resistant.
- Foot section weighs 50% less than conventional foot sections to make deliveries safe and easy.
- Once installed, the high/low shaft will stay with the bed to eliminate loss or forgetting on future installations. (Figure C)
- Labels and color-coded springs assure accurate installation of side rails for patient safety. (Figure D)
- Unique motor is completely self-contained to reduce weight and noise.
- The motor assembly can be installed and removed with patient in bed.
- Bed, side rails and mattress safety-tested by SGS, a nationally recognized testing facility. Satisfies all applicable FDA entrapment guidelines.
- The motor is UL approved.
- In the event of a power failure, a 9 volt battery, located in the motor, can lower the head and foot sections nine times. No crank needed.
- New and improved hand pendant has large, easy-to-use controls. (Figure E)
- Bed ships in two cartons.
- Channel frame construction provides superior strength and reduced weight.
- Zinc-coated spring deck.
- Extension kit available (model #15030EXTKIT) extends length from 80" to 84".
- Patent Pending.
- Weight Capacity 450 lbs.

specifications

Bed Height Minimum Without Casters	12.75"
Bed Height Maximum With Casters	24"
Unit Dimensions	36" (W) x 88" (L)
Weight	163 lbs.
Carton Shipping Weight	179 lbs.
Weight Limit	450 lbs.
Warranty	(Original Purchaser) 5 Years on Motor and Hand Control, Limited Lifetime on Welds and Frame, 1 Year on All Other Parts and Components

accessories

15030HBC	Bedding in a Box - contains 1 pillow case, 1 fitted sheet, 1 cover sheet, and 1 blanket.	1 set per bx/5 bxs per cs
15030HBL	Fitted Sheet - Contains 2 Fitted Sheets	2 sheets per bx/ 10 bxs per cs