

COMPETITIVE EDGE Line



HCPSC
E0261



15560



A



B




C




D



E

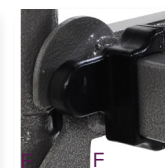

F

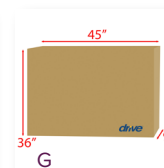

G

Competitor™ Semi-Electric Bed

Item #	Product Description	UOM
15560	Semi-Electric Bed with Full Side Rails	1/cs
15560-PKG	Semi-Electric Bed with Full Side Rails and Mattress (item #15006)	1/cs

- Receiver bracket mounts on the frame to eliminate the need for cross braces. Rail receiver brackets are retractable. Ensures proper placement and reduces lost or forgotten components. (Figure A)
- Hand pendant has large, glow-in-the-dark, easy-to-use controls. Battery backup in hand control. Also has troubleshooting light built in. (Figure B - C)
- 15560 comes in ultra-compact carton (36" x 9" x 45"); all components come in one carton. (Figure G)
- Adjustable height bed. Leg height (Figure D) and deck height (Figure E) are adjustable without tools.
- Safety lock slides over hanger pins to prevent disengaging. (Figure F)
- Motor is completely self-contained to reduce weight and noise (same as Delta™ Ultra Light 1000).
- Motor assembly can be installed and removed with patient in bed.
- Motor is UL approved.
- Bed comes complete with full side rails and retractable rail sockets to mount rails.
- Deck height of the bed can easily be adjusted for multiple heights without tools.
- Casters are an option.
- Head and foot frame weigh only 42 lbs. each.
- Head and foot boards are interchangeable.
- Optional trapeze available (model #15560TRAPEZE).

specifications

Bed Height Minimum Without Casters	11" - 20"
Bed Height Maximum With Casters	14" - 23"
Unit Dimensions	35" (W) x 84" (L)
Weight	143 lbs.
Carton Shipping Weight	151 lbs.
User Weight Limit	450 lbs.
Warranty	(Original Purchaser) 5 Years on Motor and Hand Control, Limited Lifetime on Welds and Frame, 1 Year on All Other Parts and Components

accessories

15560TRAPEZE	Trapeze use with Competitor™ Bed (15560) (2 locking, 2 non-locking)	1/cs
15044CSET	Set of 4 Casters for use with Competitor™ Bed (15560)	4/pkg