

* HCPCS CODING:
E2603 / E2604

skin protection




ASCENT

*see chart for dimensions

PELVIC, LATERAL, AND MEDIAL THIGH
SUPPORTS PROVIDE
OPTIMAL POSITIONING.

VISCO INSERT PROVIDES
EXTRA PROTECTION
FOR ISCHIAL AREAS.

TAPERED LEG ADDUCTORS
INCREASE IN WIDTH NEAR THE
BACK OF THE CUSHION TO GIVE
SUPPORT TO THE GREATER
TROCHANTER SHELVES.

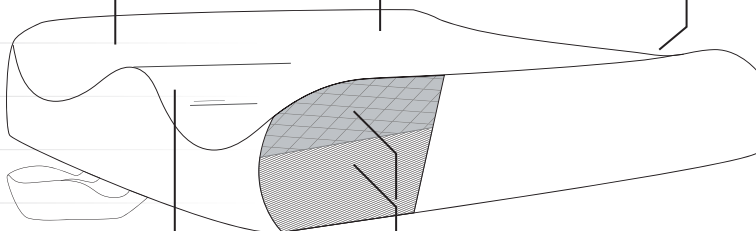


MEDIAL THIGH SEPARATOR
PROMOTES PROPER
LEG ALIGNMENT.

TWO DENSITIES OF HIGH
RESILIENCY FOAM INCREASE
SURFACE AREA AND COMFORT.

VISCO ISCHIAL INSERT

VISCO FOAM CONTOURS TO
USERS BODY, PROVIDING
MAXIMUM SKIN PROTECTION.

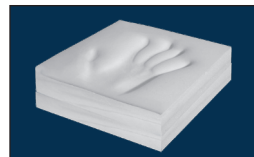

-covers-

STRETCH-AIR FABRIC — designed for
air transmission and heat reduction. This
two-ply, breathable laminate combines a high-
stretch polyester top layer with a waterproof,
breathable vapor barrier.

COMFORT-TEK — fluid proof exterior has
supple leathery properties that provide an
easily cleaned surface and protection for the
inner cushion. Comfort-Tek has multidimensional
stretch and helps alleviate pressure on bony
prominences.

CUSHION STYLE: Visco Insert		FABRIC TYPE: Stretch-Air			PRODUCT TYPE: Seating Surface		
MODEL NUMBER	W X L	FRONT	REAR	ABDUCTOR	LEG TROUGH	WT. LIMIT	
AC-GS-1616	16 x 16	3"	2.75"	3.5"	2.5"	250 lbs.	
AC-GS-1618	16 x 18	3"	2.75"	3.5"	2.5"	250 lbs.	
AC-GS-1620	16 x 20	3"	2.75"	3.5"	2.5"	250 lbs.	
AC-GS-1816	18 x 16	3"	2.75"	3.5"	2.5"	300 lbs.	
AC-GS-1818	18 x 18	3"	2.75"	3.5"	2.5"	300 lbs.	
AC-GS-1820	18 x 20	3"	2.75"	3.5"	2.5"	300 lbs.	
AC-GS-2016	20 x 16	3"	2.75"	3.5"	2.5"	350 lbs.	
AC-GS-2018	20 x 18	3"	2.75"	3.5"	2.5"	350 lbs.	
AC-GS-2020	20 x 20	3"	2.75"	3.5"	2.5"	350 lbs.	
AC-GS-2216	22 x 16	3"	2.75"	3.5"	2.5"	400 lbs.	
AC-GS-2218	22 x 18	3"	2.75"	3.5"	2.5"	400 lbs.	
AC-GS-2220	22 x 20	3"	2.75"	3.5"	2.5"	400 lbs.	
AC-GF-2416	24 x 16	3"	2.75"	3.5"	2.5"	400 lbs.	
AC-GF-2418	24 x 18	3"	2.75"	3.5"	2.5"	400 lbs.	
AC-GF-2420	24 x 20	3"	2.75"	3.5"	2.5"	400 lbs.	

ALL CUSHIONS AVAILABLE WITH RIGID INSERT AND EXTRA INCONTINENT PROOF LINER (PAGE 11)

CUSHION STYLE: Visco Insert		FABRIC TYPE: Comfort-Tek			PRODUCT TYPE: Seating Surface		
MODEL NUMBER	W X L	FRONT	REAR	ABDUCTOR	LEG TROUGH	WT. LIMIT	
AC-GF-1616	16 x 16	3"	2.75"	3.5"	2.5"	250 lbs.	
AC-GF-1618	16 x 18	3"	2.75"	3.5"	2.5"	250 lbs.	
AC-GF-1620	16 x 20	3"	2.75"	3.5"	2.5"	250 lbs.	
AC-GF-1816	18 x 16	3"	2.75"	3.5"	2.5"	300 lbs.	
AC-GF-1818	18 x 18	3"	2.75"	3.5"	2.5"	300 lbs.	
AC-GF-1820	18 x 20	3"	2.75"	3.5"	2.5"	300 lbs.	
AC-GF-2016	20 x 16	3"	2.75"	3.5"	2.5"	350 lbs.	
AC-GF-2018	20 x 18	3"	2.75"	3.5"	2.5"	350 lbs.	
AC-GF-2020	20 x 20	3"	2.75"	3.5"	2.5"	350 lbs.	
AC-GF-2216	22 x 16	3"	2.75"	3.5"	2.5"	400 lbs.	
AC-GF-2218	22 x 18	3"	2.75"	3.5"	2.5"	400 lbs.	
AC-GF-2220	22 x 20	3"	2.75"	3.5"	2.5"	400 lbs.	
AC-GF-2416	24 x 16	3"	2.75"	3.5"	2.5"	400 lbs.	
AC-GF-2418	24 x 18	3"	2.75"	3.5"	2.5"	400 lbs.	
AC-GF-2420	24 x 20	3"	2.75"	3.5"	2.5"	400 lbs.	

ALL CUSHIONS AVAILABLE WITH RIGID INSERT AND EXTRA INCONTINENT PROOF LINER (PAGE 11)